
Ulike

Air 10 IPL Hair Removal Device User Manual

Air Series

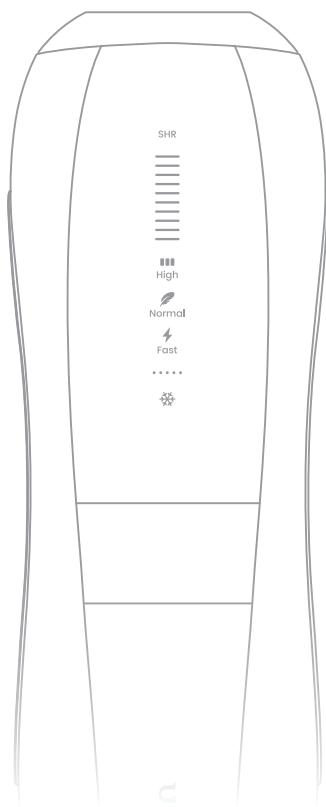ENGLISH ————— 02-19

FRANÇAIS ————— 22-39

ESPAÑOL ————— 42-59

Thank you for choosing Ulike.

Founded in 2013, Ulike has focused on at-home IPL hair removal devices from day one.

We believe hair removal should be a safe, smooth, and comfortable experience.

And this device is going to be the hair removal system you like.

Feel free to contact us via email, telephone, or chat with us through our website if you need assistance or have any questions. We're here to help!

Official Website: support@ulikeglobal.com

Channel Support: contact@ulike.com

Tel: +1 855-577-8107

Amazon Support: amazonsupport@ulike.com

Tel: +1 833-206-6933

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www.ulike.com

Shenzhen Ulike Smart Electronics Co., Ltd.

810, Building 1, Xunmei Science and Technology Plaza, No. 8 Keyuan Road, Science Park Community, Yuehai Sub-District, Nanshan District, Shenzhen 518000, Guangdong, P.R. China

Before using the device, please read the user manual carefully, including all contraindications, warnings and safety information, and keep it for future reference. Illustrations in this manual are for reference only.

Safety Precautions

■ Contraindications

Not everyone is suitable to use the IPL Hair Removal Device, please refer to the following information in detail before using.

- DO NOT use the device on skin that has burns, tattoos, or scars to prevent injury.
- Do NOT use the device if your skin has been artificially or naturally tanned recently or burned from overexposure to sunlight. Your skin may be extra sensitive following exposure and particularly susceptible to side effects of IPL treatment (e.g. skin burns, blistering, discoloration, or scarring).
- DO NOT use if you are sensitive to sunlight (photosensitivity) or are taking medication that causes photosensitivity e.g., certain antibiotics such as tetracycline, NSAIDS such as ibuprofen or retinoids such as Accutane and/or topical retinoids such as Retin A. Always check the instructions that come with your medical products to see if photosensitivity is listed as a side effect. Seek for medical advice if you are unsure of whether it is suitable to use IPL Hair Removal Device. Using the device on photosensitive skin could result in skin damage including blisters and swelling.
- DO NOT use the device if you have a history of skin cancer or pre-cancerous lesions (e.g. nevi or a large number of moles).
- DO NOT use if you are pregnant, lactating or menstruating.
- DO NOT use the device on people under the age of 18.

■ Warnings

- DO NOT use the device on the following areas:
 - Testicles, genitals, and nipples.
 - Dark brown or black spots such as birthmarks, moles or warts.
 - Tattoos or long-lasting make up.
 - Areas affected by skin conditions such as psoriasis, vitiligo, eczema, acne, herpes simplex, active infections, or wounds.
 - Areas that have had plastic surgery.
 - Near or around the eyes.
- DO NOT exceed the usage specified in this manual. Damage of the skin may occur after prolonged or repeated surface application on one site.

■ EYE SAFETY

The device emits intense pulses of light during treatment, and direct exposure can cause damage to the eyes.

- DO NOT look directly into the sapphire light window when the device is powered on.
- DO NOT attempt to activate the device towards the eyes.
- DO NOT treat areas around the eyes, such as eyebrows or eyelashes.

It is recommended to wear goggles and look away from the device when activating the lighting button. Additionally, it is advised to avoid having the air inlet facing toward your eyes during treatment. This allows for proper contact of the device with the treatment area, ensuring effective light transmission while shielding the eyes from stray light.

■ ELECTRICAL AND FIRE SAFETY

The PATIENT is an intended OPERATOR.

- DO NOT disassemble, repair and refit the product without the presence of a professional personnel as this may cause damage, malfunction or personal injury.

After-service is provided by an authorized service spot. If needed, circuit diagrams, component part lists, descriptions, calibration instructions can be provided by MANUFACTURER to assist to SERVICE PERSONNEL in parts repair.

Any unauthorized modification is prohibited, which may lead to unacceptable risk.

- DO NOT use if the device is damaged or broken; for example, cracked housing, cable damage (internal wires visible), cracked or broken light window etc.
- DO NOT use near bathtubs, showers, basins or other vessels containing water. The device should not be used if it is exposed to water.
- DO NOT use if the device becomes uncomfortably hot to touch.
- DO NOT attempt to activate (flash) the device onto any surface other than skin.
- DO NOT block the air vents on both sides of the device.
- DO NOT carry out maintenance during use.
- DO NOT leave the device unattended when it is switched on to avoid the risk of fire or burns.

■ Unsuitable Situations

- This device is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the device to avoid the risk of fire or burns.
- This device is not intended for use by children. Children should be supervised to ensure that they do not play with the device to avoid the risk of fire and burns.
- If the main cord is damaged, you should contact our support center and have it replaced by Shenzhen Ulike Smart Electronics Co., Ltd., a service center authorized by Shenzhen Ulike Smart Electronics Co., Ltd. or similarly qualified persons in order to avoid risk of electrocution.
- Water and electricity are a dangerous combination! To avoid risk of electrocution:
 - DO NOT use this device in wet surroundings (e.g. in the bathroom or near a shower or swimming pool);
 - DO NOT let water run into the appliance.
- Check if the voltage indicated on the device corresponds to the local mains voltage before you connect the device to avoid risk of electrocution or long-lasting damage to the device.
- If the wall socket used to power the device has poor connections, the plug of the device becomes hot. Make sure you plug the device into a properly installed wall socket to avoid the risk of fire and burns.
- DO NOT subject the device to heavy shocks to avoid risk of damage to the lamp.
- CAUTION: Always unplug the appliance after use and in case of a power failure to avoid risk of damage to the device.
- The risk group classification: Exempt Group.
- Keep the device away from lint, dust and light, or the performance will be affected.
- Keep the device out of children and pets' reach.
- DO NOT expose treated areas to the sun. Wait at least 48 hours after treatment before exposing treated skin to direct sunlight. Your skin may be extra sensitive following IPL treatment and particularly susceptible to sunburn. Apply sunscreen (SPF 15 or greater) to the treated area or cover the treated area with suitable clothing.

■ Precautions

- The device is intended for a single user.
- Professional skills and training are not required for the usage, please check the user manual before use.
- This device is not recommended for use on red, grey, blonde or very light hair since it is not effective on these hair colors.

Product information

Product name: **IPL Hair Removal Device**

UI20 WH, UI20 BL, UI20 PN, UI20 BR, UI20 GR, UI20 BS, UI20 MP, UI20 WG

■ Indication for Use

IPL Hair Removal Device is indicated for the removal of unwanted hair. The device is also indicated for the long-lasting reduction in hair regrowth, defined as the long-term, stable reduction in the number of hairs regrowing when measured at 6, 9 and 12 months after the completion of a treatment regime.

■ Package Contents

The package contains the following contents:

IPL Hair Removal Device*1

Adapter*1



Power Cord*1

User Manual*1



Razor*1

Goggles*1

■ Product Structure

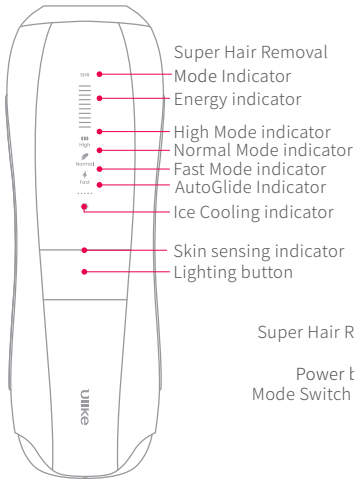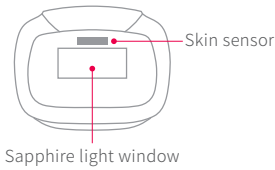IPL Device

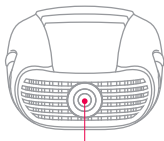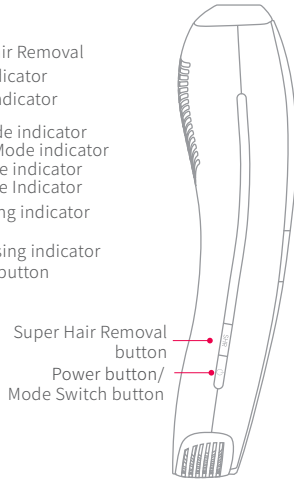Power input interface


Goggles

Adapter

Power cord

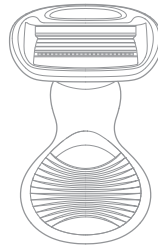Razor

■ Operation Principle

IPL Hair Removal Device is a light-based device for long-term hair removal, designed for home use. It is intended for the removal of unwanted hair and long-lasting reduction in hair regrowth. Ideal body areas include the armpits, bikini line, arms and legs. The device uses the IPL technology with lower energy level, including 4 power modes. Intense Pulsed Light technology is able to achieve long-term hair removal results at a fraction of the energy level used in other light-based hair removal devices.

The principle of operation is selective photo thermolysis and the mechanism of operation is to disable hair growth using light to preferentially heat the hair bulb.

How to Use

■ Application Ranges

The IPL Hair Removal Device is for hair removal of large areas (e.g. arms, legs, chest) and small areas (e.g. fingers, upper lip). Please check the Skin Tone Table below to determine whether IPL is suitable for you.

USE ON ACCEPTABLE SKIN TYPES AND HAIR TYPES

This device is to only be used on Fitzpatrick Skin Types I-VI. Please refer to the chart below.

- Please IDENTIFY your skin type and check your natural hair color to see whether the product is applicable:

Skin type Hair color						
	I	II	III	IV	V	VI
White/Grey	✗	✗	✗	✗	✗	✗
Red	✗	✗	✗	✗	✗	✗
Light blond	✗	✗	✗	✗	✗	✗
Light brown	✓	✓	✓	✓	✓	✗
Brown	✓	✓	✓	✓	✓	✗
Dark brown	✓	✓	✓	✓	✓	✗
Black	✓	✓	✓	✓	✓	✗

(Note: X means not applicable, ✓ means applicable; skin types are categorized according to the Fitzpatrick classification.)

NOTE:

- Select the treatment mode that best suits your skin tolerance and, if there are no adverse reactions, transition gradually to higher energy modes.
- The color of your skin may change, for example because your skin has been tanned: it may then be darker than the skin tone during the previous treatment. Therefore, it is important that you check the settings to be used before each treatment and set different settings according to

different treatment areas. Some areas may be darker than others.

- Users with skin type V should use the device with a particular caution because they are more likely to experience adverse events. If you have trouble identifying yourself as Skin Type V versus VI, please consult your physician.

■ Prior to Use

Uses the razor provided along with the device to remove all visible hair on the skin before using the device. Make sure there is no hair remaining above the surface of the skin as this can cause treatment discomfort. Remove any cosmetics, lotions and creams from the area you wish to treat. Cleanse the treatment area and pat dry.


■ Skin Test

- Before your first treatment, we recommend you test your skin for reactions.


- Select part of body you intend to treat. Connect the adapter to the device, then attach the power cord, and finally plug the power cord into the power source. Press the power button to turn it on. Press sapphire light window firmly on the selected skin until the skin sensor indicator lights up. Move the device to another test position and change treatment mode, press the lighting button to do another test treatment. If it feels comfortable, change treatment mode and try it on another test position.

- Do not apply more than one flash on the same spot. Wait 24 hours after testing, to ensure your skin is suitable for treatment and there is no adverse reaction to the light energy. If there is no abnormal reaction to the skin after 24 hours, and you didn't feel intense hotness during the treatment, then the higher treatment mode you used is suitable for your skin. You can continue to test the next energy level till you find the maximum acceptable level for the intended skin area.


■ Instructions

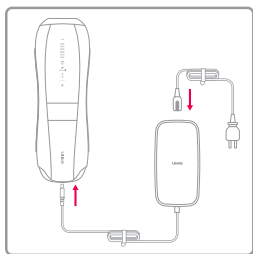
1. Connect the power cord, adapter and device as shown in figure.

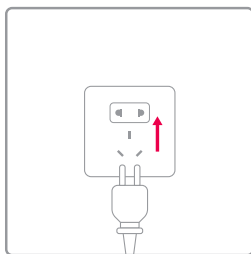2. Plug the other end of power cord into an electrical outlet. (Note: It is a plug-in device, please always keep it connected to the power during the treatment. Otherwise, it will not work).

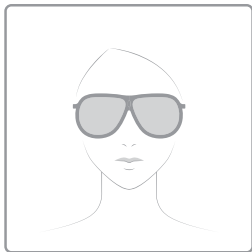
3. Wear the goggles before treatment starts.

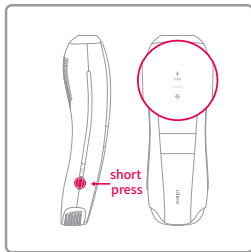4. Press the power button on the side of the device to turn it on. The default mode is Fast Mode, and the cooling function will be active during the whole process, so the relevant indicators will be on when the device is turned on.

Note: The cooling function is designed for a more comfortable hair removal experience, but the function itself doesn't have any therapeutic effect.

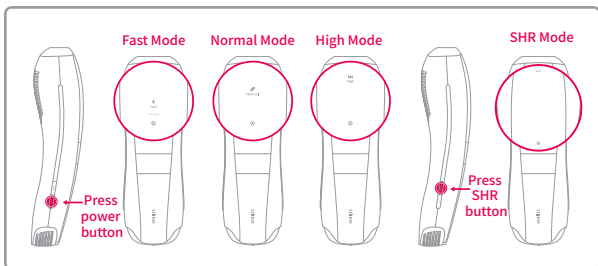5. The device offers 4 modes.

While the device is powered on, press the power button to switch to the desired treatment mode. The modes cycle in the following order: Fast Mode - Normal Mode - High Mode - Normal Mode - Fast Mode. Press SHR button to turn on/off SHR Mode.

■ Instructions for SHR Mode*

SHR Mode*: Super Hair Removal Mode

Recommend Treatment Area: Armpits and bikini area (areas with coarser hair)


Recommended Flashes on Each Treatment Area: 2-3 times

Note: SHR Mode*(*Refers to a mode that utilizes heat accumulation from multiple IPL pulses.)

Tip:


When starting to use the device, try the FAST mode for a skin test first. After confirming that there are no abnormalities in the treated area, you can use this product with confidence in Fast—Normal—High—SHR modes.

Due to the high energy of SHR, please use other three treatment modes first, and then choose the SHR mode for hair removal after your skin has established tolerance and all intended hair removal areas have been patch tested for reactivity.

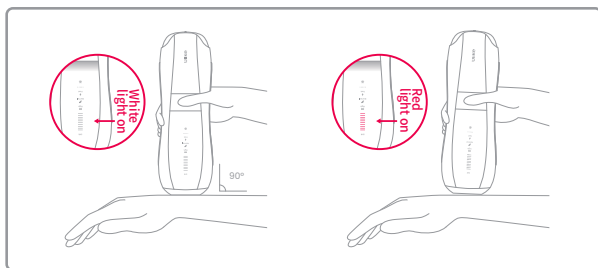6. **A,** Firstly, position the sapphire light window vertically on the area of skin to be treated,
- B,** Ensure full contact with your skin, indicated by the skin sensing indicator illuminating white.
- C,** Lastly, press the lighting button to initiate hair removal.


Note: If the red light illuminates, it indicates insufficient contact between the sapphire sapphire light window and the skin, and the device will not emit a light pulse.

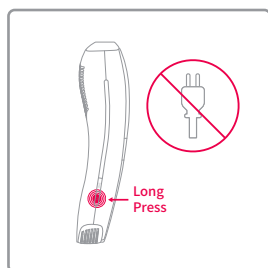7. After finishing the treatment, press and hold the power button for 1.2 seconds in order to turn off the device. The indicator light will turn off, signaling that the device is powering down.
It's important not to disconnect the device directly without properly shutting it down first.

Note: The device has a fan speed detector. When the fan speed is lower than 70% of the rated speed, the device will automatically shut down, so as to avoid the local high temperature leading to burns.

■ Treatment Methods

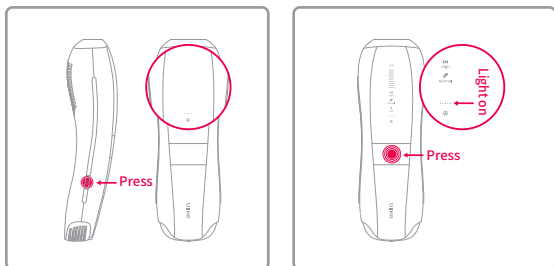1. The device offers two distinct operating modes: Auto-Glide Mode and Manual Mode.
2. How to switch the treatment methods? Long press the lighting button for 1 to 2 seconds to switch to AutoGlide Mode or Manual Mode.
3. Note, Fast Mode, Normal Mode and High Mode, support both Auto-Glide Mode and Manual Mode; SHR Mode only supports Manual Mode.


• Manual Mode (suitable for small areas such as the upper lip):

Manual mode is the default treatment method when the device powers on. Position the sapphire light window vertically on the treatment area and fully in contact with the skin. And the skin sensing indicator will light up. Then press the lighting button to start hair removal. Move the sapphire light window to the next area and repeat the steps above.

• AutoGlide (for large areas such as arms, legs):

Long press the lighting button for 1 or 2 seconds, the Auto-Glide indicator will light up. Please keep the sapphire light window in full contact with the skin and keep moving the device on the treatment area to ensure that the skin is covered evenly and avoid missing treatment areas, the device automatically emits light pulses at predefined intervals, eliminating the need for manual intervention.

Under Auto-Glide mode, press the lighting button for 1 or 2 seconds again to switch to Manual mode.

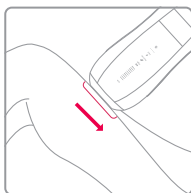■ Suggested Treatment Regime

Weeks 1-4: Three times a week, every other day.

From week 5: Adjust the treatment schedule to once every two weeks or once a month, based on your hair growth.

Please note that you can customize the treatment schedule and mode according to your hair type. For thicker or heavier hair, you can consider increasing the modes and frequency within your tolerance.

Keep in mind that individual hair growth varies, and results may differ for each person.

■ Skincare After Treatment

Following hair removal, it is advisable to apply body lotion or cold compress gel for proper skincare.

Additionally, refrain from exposing treated areas to sunlight. Wait at least 48 hours after treatment before allowing direct sunlight exposure to the treated skin. And Avoid skincare products containing alcohol and acidic/alkaline substances.

Maintenance

■ Cleaning

In order to ensure the best results and service life of the product, please clean the product before and after each use to avoid debris accumulation.

- Before cleaning, the device must be switched off and unplugged to allow the device to cool.
- Use a soft, clean and dry cloth to wipe and clean the light window surface.
- Never wash the product and accessories under faucets or in dishwashers.
- Do not use abrasive cleaners (toothpaste, etc.) or corrosive liquids (gasoline, motor oil, acetone, etc.) to clean the product.

Note: If the water does not clean well, use a cloth dipped in 75% concentration alcohol to clean the device.

When cleaning is finished, please conduct a visual inspection to check whether there is residual dirt. If the device is not visually clean, please repeat the above cleaning steps.

■ Storage

- Clean the device before you store it.
- Store the device in a dry, dust-free environment at -10°C ~+55°C / 14 °F ~ 131 °F

Disposal

In order to minimize hazards to health and the environment and ensure that materials can be recycled, this product should be disposed of at a separate collection facility for waste electrical and electronic device according to local environmental requirements.

Troubleshooting

Below summarizes the most common problems you could encounter with the device. If you are unable to solve the problem with the information below, please contact your local distributor or the seller.

Problem	Possible Reason	Solution
Unable to boot normally	Power is not connected	Check if the device is plugged in.
Unable to emit light	The sapphire light window may not be fully in contact with your skin	Ensure the sapphire light window fully covers your skin; once the skin sensing indicator illuminates, press the lighting button.
Unbearable pain during hair removal	Did not shave hair in treated area	Before using, shave the hair in the area to be treated. Failure to do so may result in the device burning surface hair, potentially causing discomfort.
	Selected treatment mode is too powerful	Select a lower-intensity treatment mode
	The light window is damaged	After 30 seconds of powering on, inspect the sapphire light window for any water mist. If none is present, please reach out to Ulike Support for further assistance.
After using, the hair removal area turns red	A slight reddening is normal and harmless to the body, and it will soon fade	Apply aloe vera gel to soothe the skin, or utilize the cooling effect of the sapphire light window.
Unsatisfactory hair removal effect	This device doesn't work on your skin tone or hair color	If your hair is pale blonde, gray, red or white, the treatment will not be effective.
	Missed some spots in the treatment area	Flashes may have been fired in close in close proximity and there should be some overlap between the flash areas.
	You use an energy level that is too low for you	Please choose a higher level next time
	Insufficient frequency of use	For better hair removal results, we recommend the following treatment cycle: Weeks 1-4: Every other day Weeks 5-12: Once every 2 weeks or once a month
	Your response to IPL light exposure is slower than expected, or the hair you'd like to remove is very thick	Thicker hair areas, such as the bikini line, armpits, and beard, may require more time and repeated treatments to see significant reduction. Different skin types may take longer to respond to light, so continued use of our product for up to 3 months is recommended for optimal results.
The device becomes hot and humid during use	May be due to prolonged operation	Turn off the device and let it stand for several minutes

Technical Specifications

Product Name	IPL Hair Removal Device
Product Model	UI20 WH, UI20 BL, UI20 PN, UI20 BR, UI20 GR, UI20 BS, UI20 MP, UI20 WG
Dimensions (W x H x D)	206.73mm*68.68mm*54.29mm
Power Input	24V $\approx$ 3.75A, 90W
Energy Level	1~10 level
Spot Size	3.9cm ² /0.6 sq in
Wavelength	550nm~1200nm
Output Energy	7.5J~26.0J
Energy Density	1.92J/cm ² ~6.67J/cm ²
Pulse Interval	0.40~3s
Software Version:	V1.0
Operating conditions	Operating temperature: 5°C ~ 30°C/ 41°F ~ 86 °F Humidity: 5%~ 80%RH Relative pressure: 80kPa ~106kPa
Storage and transportation conditions	Storage temperature:-10°C ~ 55°C/ 14 °F ~ 131 °F Humidity: 5%~ 90%RH Relative pressure: 50kPa~106kPa
Adapter	Manufacturer: ASIAN POWER DEVICES INC. Model: DA-90Q24 INPUT: 100V-240V~, 50/60Hz, 1.5A OUTPUT: 24.0V $\approx$ 3.75A, 90W
	Manufacturer: Mass Power Electronics Inc. Model: E090-1E240375M2 INPUT: 100V-240V~, 50/60Hz, 1.5A OUTPUT: 24.0V $\approx$ 3.75A, 90W
After removing the device from storage at the lowest ambient temperature (-10°C), it takes a certain amount of time for the device to return to normal operation at an ambient temperature of 20°C.	2 hours
After removing the device from storage at the highest ambient temperature (60°C), it takes a certain amount of time for the device to return to normal operation at an ambient temperature of 20°C.	2 hours

EMC & Biocompatibility

This device is Class II Medical Device. It complies with Medical Electrical Safety Standards (IEC60601-1).

This device is also complied to Medical EMC Standard (IEC 60601-1-2). All the user directly contracting materials for main device housing and output contacts in this device are biocompatible for its intended use. They are complied with biocompatibility standards ISO 10993-5 (Cytotoxicity), ISO 10993-10 (Sensitization) and ISO 10993-23 (Irritation).

EMC guidance

- **Warning:** Don't use near active HF surgical device and the RF shielded room of an ME system for magnetic resonance imaging, where the intensity of EM disturbances is high.
- **Warning:** Use of this device adjacent to or stacked with other device should be avoided because it could result in improper operation. If such use is necessary, this device and the other device should be observed to verify that they are operating normally.
- **Warning:** Use of accessories, transducers, and cables other than those specified or provided by the manufacturer of this device could result in increased electromagnetic emissions or decreased electromagnetic immunity of this device and result in improper operation.
- **Warning:** Portable RF communications device (including peripherals such as antenna cables and external antennas) should be used no closer than 12 in (30 cm) to any part of the device, including cables specified by the manufacturer. Otherwise, degradation of the performance of this device could result.

Symbols

No.	Symbols	Description
1		Indoor Use
2		"WEEE (Waste Electrical and Electronic device)". The waste products should be handled legally.
3		Keep away from rain
4		Manufacturer
5		Warning symbol for use with IPL device: Do not use it directly on the eyes, as the intense light emitted by this appliance may cause eye damage.
6		Consult instructions for use
7		Type BF Applied Part
8		DC pocket
9		Maximum number of identical transport packages/ items which may be stacked on the bottom package, where "5" is the limiting number.
10		Caution

Warranty

After-sales protection policy fully complies with relevant sales laws and regulations.

From the date of purchase, the product is covered by a 90-Day Return and 2-Year Warranty. However, warranty may be refused under the following circumstances:

- Failure to follow the instructions, operation, maintenance.
- Mistake which causes damage to the product.
- Exceeding the allowable use environment (such as high temperature, corrosion, etc.) causing product damage.
- Non-normal external forces (falling, extrusion, impact, etc.) causing damage to the product.
- Natural force majeure (fire, flood, lightning, impact, etc.) causing damage to the product.
- Use of non-manufacturers accessories and expansion parts causing damage to product.

Contact information

Shenzhen Ulike Smart Electronics Co., Ltd.
810, Building 1, Xunmei Science and Technology Plaza, No. 8
Keyuan Road, Science Park Community, Yuehai Sub-District,
Nanshan District, Shenzhen 518000, Guangdong, P.R. China

Feel free to contact us via email, telephone, or chat with us through our website if you need assistance or have any questions. We're here to help!

Official Website: www.ulike.com & support@ulikeglobal.com

Channel Support: contact@ulike.com

Amazon Support: amazonsupport@ulike.com

Appendix A: EMC Information

Essential performance:

The IPL Hair Removal Device does not have essential performance.

Guidance and Manufacturer's Declaration – Electromagnetic Emissions	
Emissions test	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class B
Harmonic emissions IEC 61000-3-2	Class A
Voltage fluctuations/- flicker emissions IEC 61000-3-3	Complies

Guidance and Manufacturer's Declaration – Electromagnetic Immunity		
Immunity Test	IEC 60601-1-2 Test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4kV, ±8 kV, ±15 kV air
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/output lines	±2 kV for power supply lines ±1 kV for input/output lines
Surge IEC61000-4-5	± 1 kV line(s) to line(s) ± 2 kV line(s) to earth	± 1 kV line(s) to line(s) ± 2 kV line(s) to earth
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % UT; 0.5 cycle at 0°, 45°, 90°, 135°, 180°, 225°, 270°, 315° 0 % UT; 1 cycle 70 % UT; 25/30 cycle 0% UT; 250/300 cycle	0 % UT; 0.5 cycle at 0°, 45°, 90°, 135°, 180°, 225°, 270°, 315° 0 % UT; 1 cycle 70 % UT; 25/30 cycle 0% UT; 250/300 cycle
Power frequency magnetic field IEC 61000-4-8	30 A/m 50Hz/60Hz	30 A/m 50Hz/60Hz
Conducted RF IEC61000-4-6	3 V rms. 150 kHz to 80 MHz 6 V RMS in the ISM and amateur bands between 0.15 MHz and 80 MHz	3 V rms. 150 kHz to 80 MHz 6 V RMS in the ISM and amateur bands between 0.15 MHz and 80 MHz
Radiated RF IEC61000-4-3	10 V/m 80 MHz – 2,7 GHz 80 % AM at 1 kHz	10 V/m 80 MHz – 2,7 GHz 80 % AM at 1 kHz

NOTE UT is the a.c. mains voltage prior to application of the test level.

Guidance and Manufacturer's Declaration - IMMUNITY to proximity fields from RF wireless communications device					
Immunity Test	IEC60601 test level				Compliance level
	Test frequency	Modulation	Maximum power	Immunity level	
Radiated RF IEC 61000-4-3	385 MHz	**Pulse Modulation: 18Hz	1.8 W	27 V/m	27 V/m
	450 MHz	*FM+ 5Hz deviation: 1kHz sine	2 W	28 V/m	28 V/m
	710 MHz 745 MHz 780 MHz	**Pulse Modulation: 18Hz	0.2 W	9 V/m	9 V/m
	810 MHz 870 MHz 930 MHz	**Pulse Modulation: 18Hz	2 W	28 V/m	28 V/m
	1720 MHz 1845 MHz 1970 MHz	**Pulse Modulation: 217Hz	2 W	28 V/m	28 V/m
	2450 MHz	**Pulse Modulation: 217Hz	2 W	28 V/m	28 V/m
	5240 MHz 5500 MHz 5785 MHz	**Pulse Modulation: 217Hz	0.2 W	9 V/m	9 V/m

Note* - As an alternative to FM modulation, 50 % pulse modulation at 18 Hz may be used because while it does not represent actual modulation, it would be worst case.

Note** - The carrier shall be modulated using a 50 % duty cycle square wave signal.

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Support canal : contact@ulike.com

Tél. : +1 855-577-8107

Support Amazon : amazonsupport@ulike.com

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